Rotomaxx[™]

High Torque, High Flow, Right Angle Mixers



Advanced Mixing Jechnologies

 $\hat{\mathbf{O}}$

Built to sanitary standards at a competitive price

For larger tank volumes requiring rapid turnover & high efficiency agitation, but also designed to meet stringent CIP washdown, pressure washing and steam cleaning conditions. The Rotomaxx is available in both a Pharmaceutical or a Food grade, all stainless configuration, with multiple output speeds and seal combinations. Our Jetfoil[™] impeller technology provides 200% of the pumping rate found on conventional axial blade impellers..

Big performance in a compact package

- Total washdown capability with (#73-01) compliance
- ✓ Tank Volumes to 10,000 gallons
- Pumping Rates to 20,000 GPM
- 1/2 to 3 HP, 30-235 RPM output speeds
- Right angle design for low headroom space
- Available with direct readout of shaft speed through internal or external proximity sensors
- ✓ Validation documentation & factory acceptance testing (FAT) for GMP and FDA requirements

With a small footprint and low headroom design, the Rotomaxx is ideal where space is critical on large tanks



Your Mixing Technology Partner . . . www.admix.com

Rotomaxx Specifications


The Rotomaxx series of stainless steel mixers are available in three (3) basic configurations as shown below:

| Mixer Model | 125 RMAX | 150 RMAX | 200 RMAX |
|----------------------------|--------------|--------------|--------------|
| Output Speed (RPM) | 30-235 | 30-235 | 30-235 |
| HP Range | 1/2 - 2 | 1/2 - 3 | 3/4 - 3 |
| Shaft Diameter | 1-1/4" | 1-1/2" | 2" |
| Maximum Shaft Length | up to 95" | up to 106" | up to 165" |
| Foil Diameter | up to 55" | up to 65" | up to 67" |
| Pumping Rate @ 1 cps (GPM) | up to 13,000 | up to 20,500 | up to 23,800 |

Rotomaxx H 125RMAX 8.7 1.25 18 .00 3.3 9.4 150RMAX 9.1 1.50 19 4.4 12 .25 200RMAX 1.00 10.4 2.00 20 5.5 12 11.8 70

Dimensions are in inches

C & E may vary with motor selection

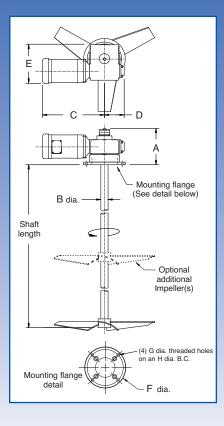

For more information:

EMPLOYEE OWNED... CUSTOMER FOCUSED

Ease of Installation & Service

The modular design of the Rotomaxx components, along with the hollow quill bore of the reducer, allows for easy mounting and servicing.

SS Motors are directly coupled to eliminate misalignment.

SS Gear Reducers include a hollow bore to eliminate shaft couplings.

SS Pedestals ..

allow easy seal disassembly and are registered for perfect seal alignment.

Removable SS Shaft

is both keyed and retained above the hollow quill to eliminate slip.

SS Jetfoil™ Impellers are keyed, welded or can


are keyed, welded or can be bolted construction where space is limited.

ADMIX

Selection & Performance

Selection of the proper Rotomaxx model for your application is simplified with our **MixxScale** sizing criteria. Our **MixxScale** program utilizes bulk fluid velocity and impeller size (foil diameter to tank diameter ratio) to provide the desired level of agitation, suspension of solids and top-to-bottom flow control. Our high efficiency Jetfoil[™] impellers provide maximum performance using the lowest horsepower consumption and an optimized gearbox size to reduce overall costs.

Our **MixxScale** (MS) selections range from 1 through 4 with MS = 1 providing a lower level of agitation with little vortex and longer blend times (10-30 minutes), while MS = 4 is characterized by very vigorous agitation with excellent suspension and very rapid blend times of under 10 minutes. The following tables provide the best Rotomaxx selection for the designated level of agitation, tank volume and viscosity.

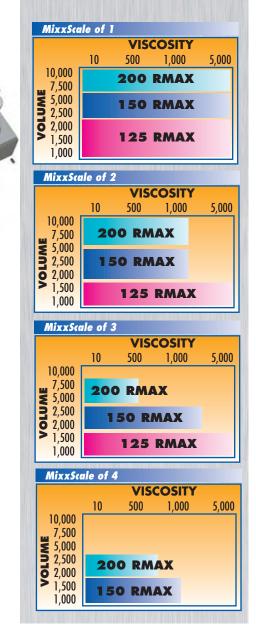

Printed in USA ROTOMAXX-2/15